




分步释放 安全可控

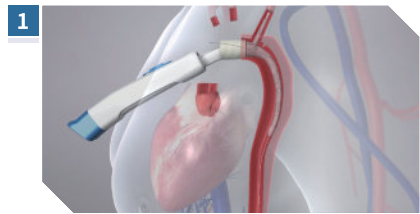


1 导入Fontus®
支架系统

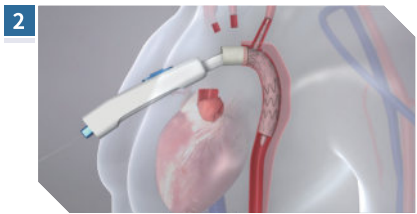

2 释放主体支架

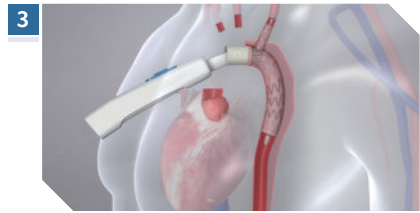


3 释放分支支架

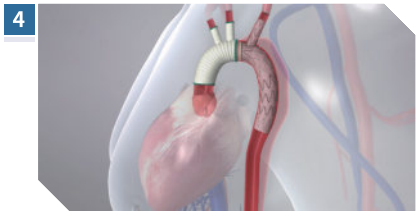

4 支架吻合

尺寸规格

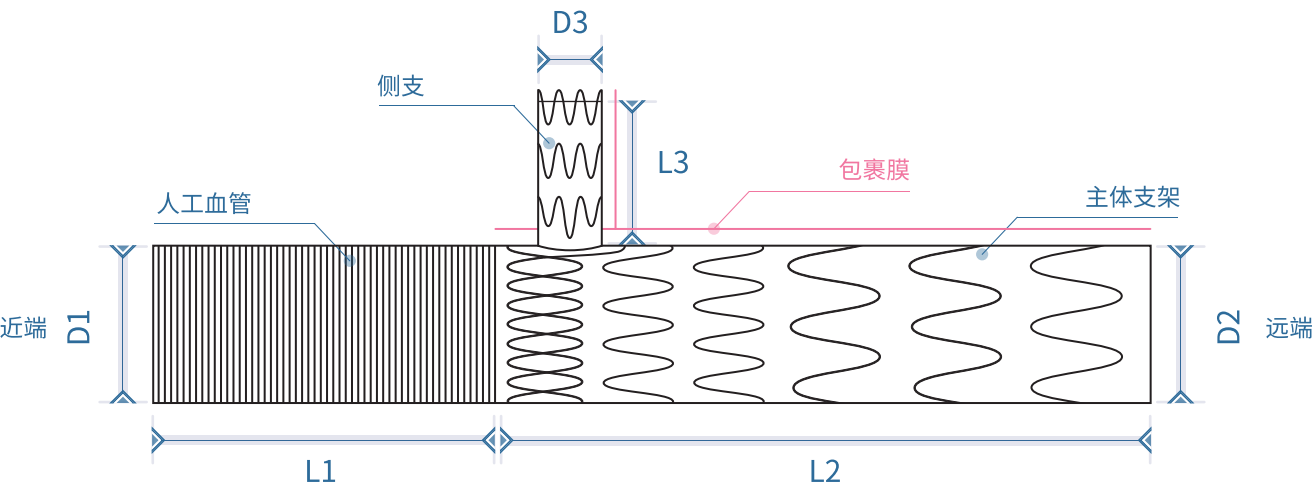

	人工血管内径 D1(mm)	主体支架远端内径 D2(mm)	侧支远端内径 D3(mm)	人工血管长度 L1(mm)	主体支架长度 L2(mm)	侧支长度 L3(mm)
FND1D2D3-L1L2L3 人工血管无涂层型号	20 22 24 26 28 30	18 20 22 24 26 28 30 32	6 8 10 12 14	10 30 60 100	80 100 120 140 160 180 200	20 25 30 35 40
FD1D2D3-L1L2L3 人工血管带涂层型号	32 34					

注:0≤D1-D2≤2




地址：上海市·浦东新区·康新公路3399弄·1号楼（上海国际医学园区·医创园）
网址：www.endovastec.com
电话：(86) (21) 38139300 传真：(86) (21) 33750026 邮编：201318




Fontus®
分支型术中支架系统

象由心生
再续经典


Fontus[®] —— 新一代术中支架系统

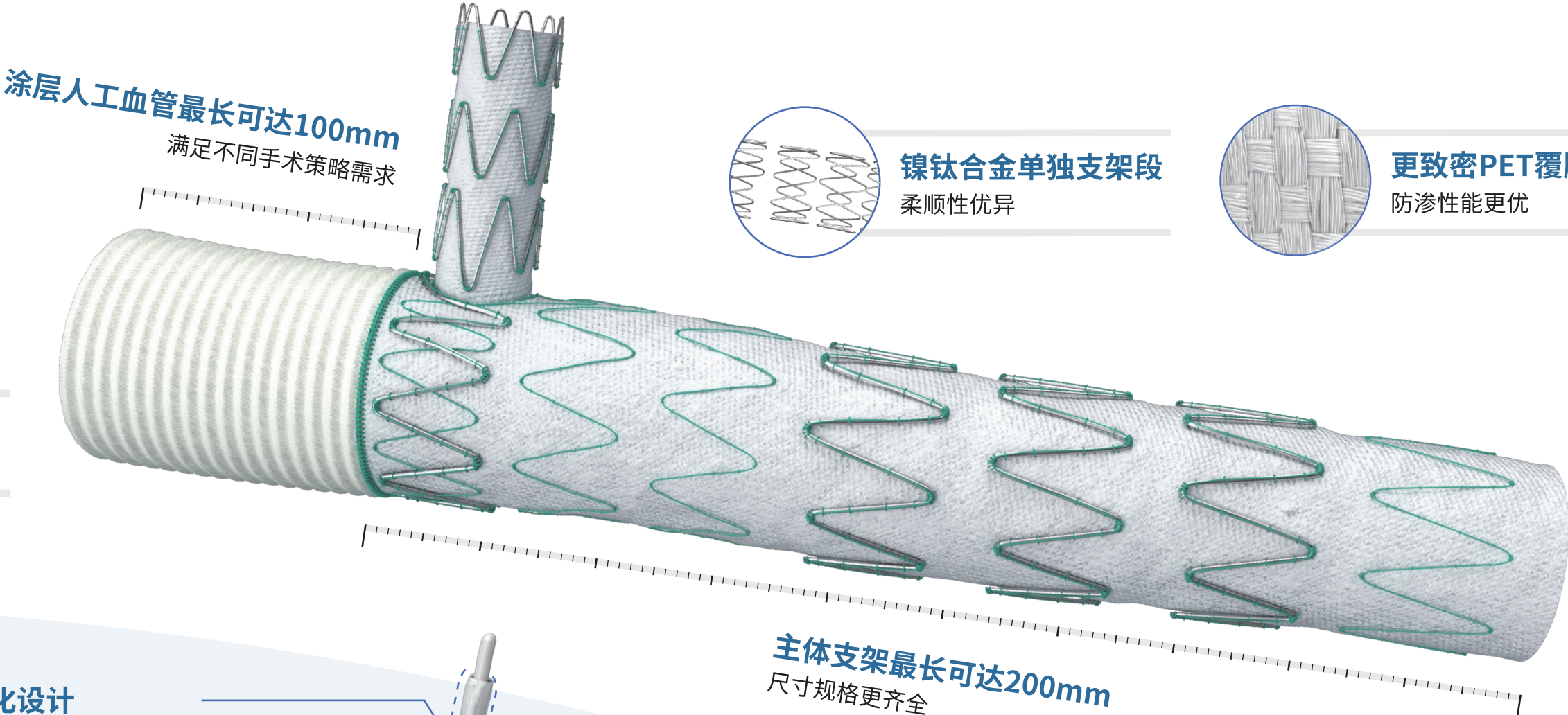


涂层人工血管最长可达100mm
满足不同手术策略需求

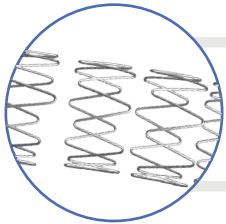

镍钛合金单独支架段
柔顺性优异

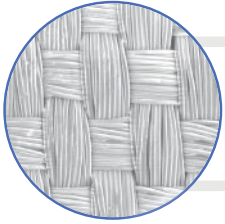


更致密PET覆膜材料
防渗性能更优

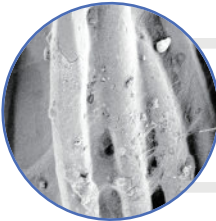


胶原蛋白涂层
防止渗血

主体支架最长可达200mm
尺寸规格更齐全

侧支一体化设计
主体吻合口前移,降低手术难度

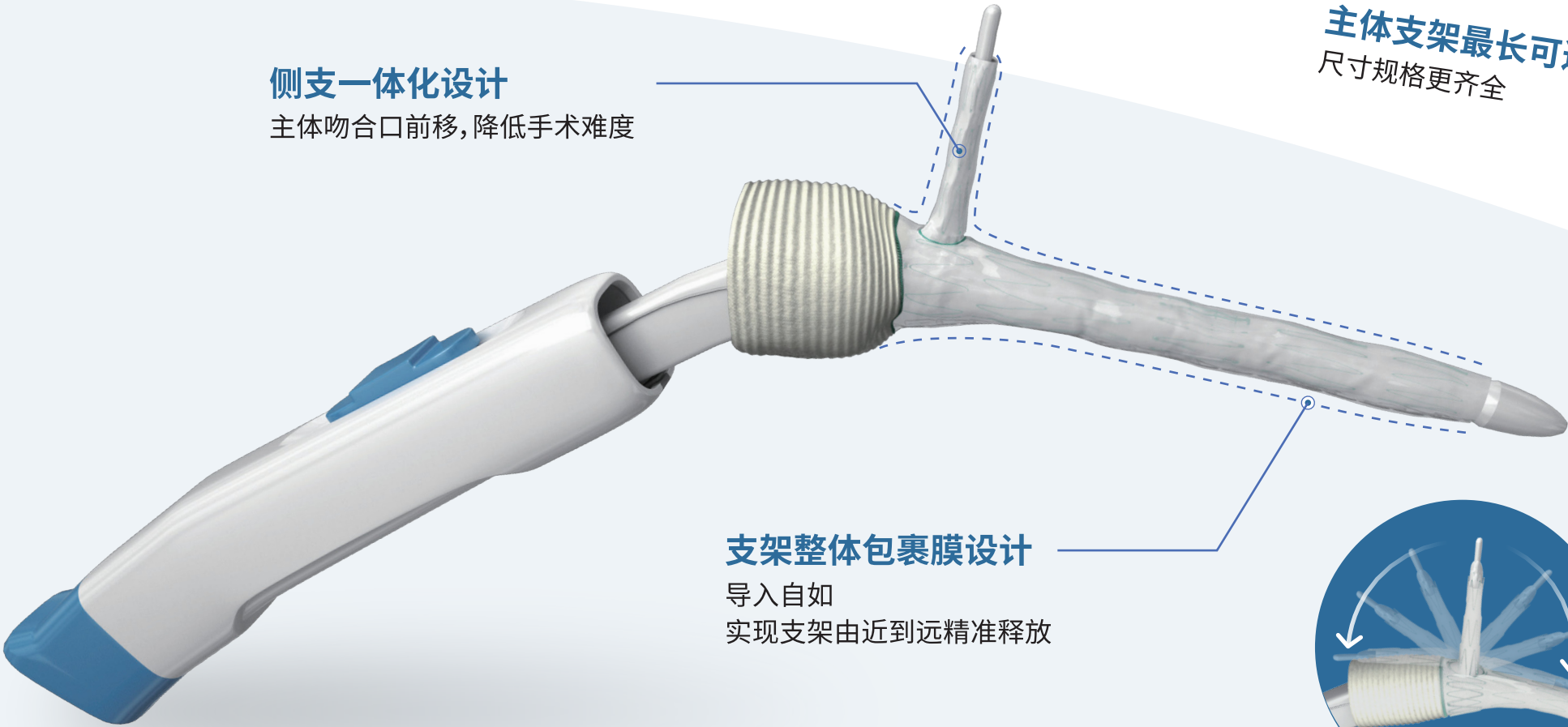

支架整体包裹膜设计
导入自如
实现支架由近到远精准释放

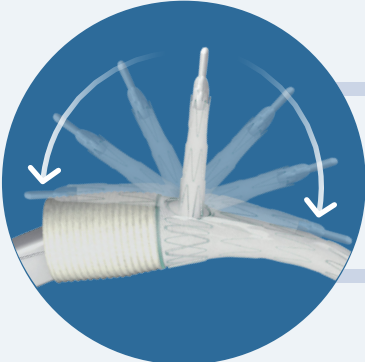

侧支角度可调
合理顺应血管方向

上一代

CRONUS[®]
术中支架系统

